

clarisonic

Opal SONIC
SKIN
INFUSION

User Guide

Sonic Infusion Technology for
Anti-Aging

Anti-Aging Sea Serum

Universal Voltage Charger

2 Replacement Applicator Tips

Two-Year Warranty

Infuses serum more effectively
than manual application for
younger looking eyes.



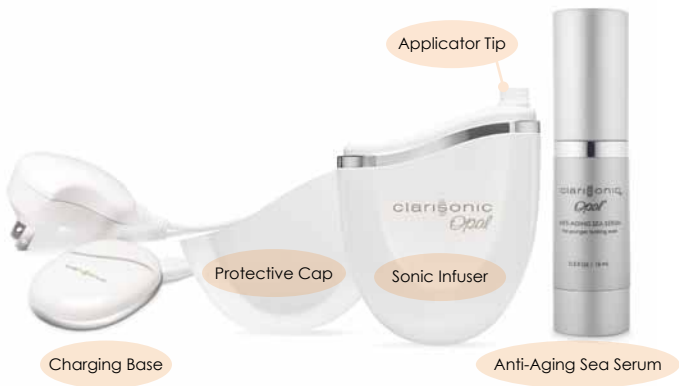


Figure 1

OPAL® SONIC INFUSION SYSTEM

The rechargeable Opal Sonic Infusion System consists of the Sonic Infuser with Applicator Tip, Anti-Aging Sea Serum, a Protective Cap, and a Charging Base.

The Sonic Infuser uses a gentle micro-massaging action to effectively, but gently, apply Serum to the area(s) you need it most. The sonic action helps infuse the Serum into the outer layers of the epidermis, where it starts working immediately to reduce the appearance of fine lines and puffiness. With daily use, the Opal Sonic Infuser can provide visibly firmer, more hydrated, younger-looking skin.

THE VALUE OF SONIC INFUSION

It is believed that anti-wrinkle lotions are more effectively integrated into the skin when applied with finger tapping. The Sonic Infuser substitutes state-of-the-art sonic technology for finger-tapping to ensure effective, consistent infusion of anti-aging treatment products into the layers of the skin where fine lines reside. The Sonic Infuser's gentle kneading action allows the Serum to be infused into the outer layers of the epidermis (shown as blue area in Figure 2). As a result, the Serum can work immediately and more effectively to help promote smooth, healthy, hydrated skin.*

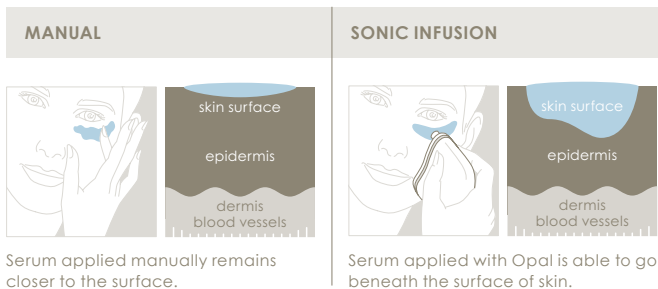


Figure 2

*Data on file. Tested using the Anti-aging Sea Serum and proven to penetrate the outer layers of the epidermis.

ANTI-AGING SEA SERUM

Robb Akridge, PhD, developed this proprietary formula with marine extracts and botanicals that help to:

Hydrate – Natural moisturizing factors and Trehalose from desert plants increase moisture content to restore natural skin hydration

Firm and tighten – Glycosaminoglycans and Cyclopeptide-5 help promote noticeably firmer and tighter skin

Protect – the delicate skin near the eye and on the face with multiple ingredients, including Canadian Willowherb and Hydroxyresveratrol

Benefits include:

- Leaves skin feeling energized, more radiant and refreshed
- Immediate and long-lasting hydration
- Skin feels firmer and tighter
- Reduces the appearance of fine lines and under-eye puffiness



DANGER

READ ALL INSTRUCTIONS BEFORE USE

IMPORTANT SAFEGUARDS

When using electrical products, basic safety precautions should always be followed, especially when children are present, including but not limited to the following:

WARNING – to reduce the risk of burns, electrocution, fire or injury to persons

- Close supervision is necessary when this product is used by, on, or near children or physically challenged persons.
- Do not use while sleeping or drowsy.
- Do not use Charging Base outdoors.
- Do not use Charging Base where oxygen is being administered.
- Keep cord away from heated surfaces.
- Do not operate Charging Base if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Call Clarisonic at 888-525-2747.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

WHEN DO I USE OPAL?

After cleansing, treat your eye with the Opal.

Opal is gentle enough for daily use, morning and evening. It's designed to target fine lines and puffiness.

OPAL SEAMLESSLY INTEGRATES INTO YOUR ANTI-AGING RITUAL

Step 1

Cleanse

Step 2

Infuse serum with Opal instead of manual application

Step 3

If desired, layer with your favorite treatment products

IMPORTANT CAUTIONS AND WARNINGS

Consult physician before using if:

- your eyes or skin are sensitive to skin care products,
- if you think you may have a skin condition, or
- if you are under a doctor's care for a skin condition.

Refer to our website for more information.

This product uses magnets.

Do not share Applicator Tips with others.

Do not use the Opal directly on the eye or the eyelid.

Do not use on open sores or pustular acne.

PRE-USE INSTRUCTIONS

BEFORE USING CLARISONIC OPAL FOR THE FIRST TIME, PLEASE DO THE FOLLOWING:

1. ALLERGY TEST

As with any skin care product, test for skin allergy at least 48 hours before using

Apply a drop of Serum to the inside of the forearm.

Rub excess Serum into arm.

If bumps, burning, skin discoloration or other allergic reaction occurs, discontinue use.

Consult your physician if symptoms persist.

2. CHARGING THE OPAL SONIC INFUSER

Charge the Sonic Infuser for 24 hours. A fully charged Sonic Infuser provides 21 days of use when used twice per day, for two 30-second uses. To safely charge or recharge the Sonic Infuser, read all warnings and instructions below.

DANGER – Electrocutation Risk With Charger

- Keep charger away from water.
- **If charger falls in water, do not touch it!** Unplug from wall immediately.
- Never use while bathing.
- Do not use or store where charger can fall or be pulled into tub or sink.

Charging/Recharging Instructions

Make sure the Protective Cap is attached to the Sonic Infuser.

Attach the power cord to the charger (factory setting has power cord attached).

Plug the cord into a working outlet.

Place the Opal Sonic Infuser on the Charging Base. Align the stripes of the Charging Base with the stripes of the Sonic Infuser, with Clarisonic logos facing up as shown in Figure 3. Feel a magnetic click.

When the Opal Sonic Infuser is properly placed on the Charging Base, it will pulse twice. A light ([the LED indicator](#)) will flash to indicate that the unit is charging. If you do not see the LED indicator, it is not properly attached.

Charge for 24 hours.



Figure 3

Charge Indicators

Action	Indicator
Attach the Opal Sonic Infuser to the charging base.	Opal Sonic Infuser pulses twice.
Charge the Sonic Infuser.	LED indicator flashes while charging. LED indicator glows constant when fully charged.
After about 40 minutes of use, the Sonic Infuser begins to indicate charge level is low.	Infuser lightly pulses and LED indicator flashes after the run cycle.
After 45 minutes of use, the Sonic Infuser needs to be charged.	Does not power on. Place on the charging base and allow 24 hours for a full charge.

The charger is compatible with 100-240 VAC; 50/60 Hz, allowing it to be used worldwide. If the shape of the plug does not fit the power outlet, use an attachment plug adapter of the proper configuration for the power outlet.

INFREQUENT USE: the Clarisonic Opal should receive a charge occasionally. If you rarely use the Opal, charge it every 3 months to maintain its charging capability.

USING THE OPAL SONIC INFUSER



Figure 4

1. APPLY SERUM TO APPLICATOR TIP

- **Separate** Cap from Sonic Infuser



- **Position the Anti-Aging Sea Serum bottle over the Applicator Tip** of the Sonic Infuser.



- **Fill** the Applicator Tip with serum.





2. APPLY SERUM

Use the Sonic Infuser around the eye area and lip contour for more effective application of serum



- **Press and release the Power button** on the side of the Sonic Infuser to activate
- **Place the activated Sonic Infuser** on the crow's feet area
- Gently **follow the rim of bone** outlining the eye area (tip should glide easily over the skin)
- Use small **overlapping circular** motion (automatically turns off after 30 seconds)
- **DO NOT USE** on soft tissue under eye and eye lids to avoid eye injury
- Gently **tap** any **remaining Serum** into the skin using the ring finger
- **Pump additional Serum** onto the Applicator Tip and **repeat for other eye**

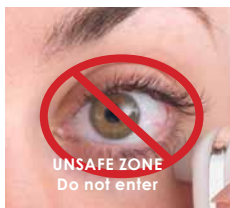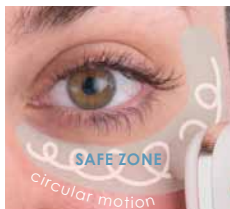

*As with any skincare product, check for skin allergies by placing a small amount to inside of forearm. If bumps, burning, skin discoloration or other allergic reaction occurs, discontinue use.

Apply Serum only to the **Safe** zone – Figure 5 shows where it is safe to use the Sonic Infuser.



Figure 5

! WARNING

To prevent eye injury, keep Applicator Tip away from eyeball and eyelids.

Use only in **SAFE ZONE**, shown above.

NOTES

- Pressing too hard will stop the motion of the Applicator Tip.
- If you press the Power button while Opal is in motion, the Infuser will turn off. To start another 30-second cycle, press the Power button.

CLEANING THE OPAL SONIC INFUSER

1. After use gently blot the Applicator Tip with a tissue or cotton swab.

DO NOT:

- Put any part of the Opal (including Applicator Tip) in water or dishwasher. **The Opal is not waterproof.**
- Use abrasive cleaners or chemicals on the Opal Sonic Infuser.
- Attempt to open the Opal Sonic Infuser. It is a sealed unit.

2. Replace Applicator Tip every 3 months.

Two replacement tips are included with this System. Tips are also sold separately on clarisonic.com, at various retail locations and from your skin care professional.

Remove Old Applicator Tip

- Using fingernail or thumb, pop off old Applicator Tip by gently grasping below the hard plastic rim. (See Figure 6.)
- To avoid damaging the tip, be careful not to pull on the soft silicone. If silicone separates from the plastic base, do not use damaged tip.

WARNING

The Applicator Tip is a small part that can cause choking.

Keep Applicator Tips away from small children.



Figure 6

Attach New Applicator Tip

- Remove new Tip from clear plastic Tip holder by grasping “wings” of plastic holder and gently pulling to break side seals. Tip will release from holder. Dispose of or recycle holder.
- Firmly press new Applicator Tip onto black tip base.
- Listen for it to snap into place to ensure proper attachment.

STORING THE OPAL

The Clarisonic Opal can be safely stored between 32° and 120°F. Subjecting Opal to extreme heat or cold could damage the unit. When not in use, it is safe to store the Opal on the Charging Base. You may keep the Charging Base plugged in if you choose, or you may disengage the power cord from the Charging Base and rest the Opal on the base. Re-engage the power cord when the Opal is ready to be recharged.

SAVE THESE INSTRUCTIONS

The Opal Sonic Infuser, Protective Cap and Charging Base have no serviceable parts.

The supply cord cannot be replaced. If the cord is damaged, do not use. Call Clarisonic with questions at 888-525-2747.

NOTE: This equipment has been tested and found to comply with Part 18 of the FCC Rules. Refer to our website for further information.

ANTI-AGING SEA SERUM INFORMATION

INGREDIENTS: Water (Aqua), Glycerin, Butylene Glycol, Chondrus Crispus (Carrageenan) Extract, Epilobium Angustifolium Flower/Leaf/Stem Extract, Inula Crithmoide Flower/Leaf Extract, Laminaria Digitata Extract, Pichia/Resveratrol Ferment Extract, Laurdimonium Hydroxypropyl Hydrolyzed Soy Protein, Trehalose, Panthenol, Sodium Hyaluronate, Glycosaminoglycans, Ascorbyl Glucoside, Ergothioneine, Disodium EDTA, Lecithin, Ferulic Acid, Hydroxyresveratrol, Hydroxyphenyl Propamidobenzoic Acid, Cyclotetrapeptide-24 Aminocyclohexane Carboxylate, Acrylates/C10-30 Alkyl Acrylate Crosspolymer, Glyceryl Acrylate/Acrylic Acid Copolymer, Alcohol, Ectoin, Sodium Hydroxide, Ethoxydiglycol, Pentylene Glycol, Phenoxyethanol, Caprylyl Glycol, Ethylhexylglycerin, Hexylene Glycol.

Anti-Aging Sea Serum Bottle 0.5 fl Oz | 15 mL

Gentle for daily use. Paraben free. No animal testing.

Serum formulated in the USA. Dermatologist and Ophthalmologist tested. Data on file.

FEEL AND SEE GUARANTEE

We are happy to offer a 90-day money back guarantee. If you are for any reason dissatisfied with the purchase of your Clarisonic Opal from an authorized dealer, call Customer Service at **888-525-2747** to make arrangements for the return.

In addition, every Clarisonic Opal comes with a two-year warranty covering any problems resulting from everyday use of the product. Call Customer Service at **888-525-2747** if you feel your Opal Sonic Infusion System is not operating properly.

For more information regarding the Opal Sonic Infusion System, visit **clarisonic.com/opal/support**.

PRODUCT REGISTRATION

Register your Clarisonic Opal so we can best serve you with product support. Plus, you can sign up for Clarisonic Email and gain access to sneak previews, new product introductions, and great insider skin care tips. Visit **clarisonic.com/registration** to register.

For more comprehensive information, new products, product user guide updates and to stay connected, visit **clarisonic.com**.

CLARISONIC® CUSTOMER SERVICE AND LIMITED WARRANTY INFORMATION

Congratulations on your purchase! We value our customers and aim to provide the highest quality products to ensure customer satisfaction. Please use this as your reference for any customer service or warranty needs you may have in the future. We also encourage you to register your product at **www.clarisonic.com/registration**.

We are happy to offer a 90 day money back guarantee. If you are for any reason dissatisfied with the purchase of your Opal from an authorized dealer, call Customer Service at 888-525-2747 to make arrangements for the return.

Clarisonic's warranty obligations are limited to the terms set out below.

1. Warranty Coverage

- a) For the duration of the Warranty Period, Clarisonic warrants the Product against defects in materials and workmanship arising from Normal Use of the Product.
- b) If there is a defect in the Product which is covered by this Limited Warranty, Clarisonic agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of warranty.
- c) This Limited Warranty is exclusive of all other warranties, whether oral or written, express or implied. There are no implied warranties created by the manufacture, sale, or use of the Product.

If Clarisonic cannot lawfully disclaim implied warranties under this Limited Warranty, all such warranties, including warranty of merchantability and fitness for a particular purpose, are limited in duration of this Limited Warranty. Some States do not allow limitations on how long an implied warranty lasts, so this limitation may not apply to you.

2. How to get warranty service or for general customer service inquiries

In order to receive service under this Limited Warranty, please call Customer Service at 1-888-525-2747 who will direct you how to proceed.

3. Limitations

- a) This warranty is limited and may be relied upon only:
 - (i) by the original end user of Clarisonic's Product;
and
 - (ii) where the Product was manufactured by or for Clarisonic and sold by an authorized distributor.
- b) To obtain service under this Limited Warranty, you will be required to provide Clarisonic with original proof of purchase date or you must have registered your Product. Registering your Product permits Clarisonic to contact you in the event of a product question or recall, allows you to sign up to receive future correspondence from us and gives you easy access to the latest Clarisonic news and information.
- c) This Limited Warranty does not apply where the Product is subjected to use that does not constitute Normal Use, including but not limited to personal injury or property damage arising from misuse of the Product.

4. Exclusions

Clarisonic is not responsible for any indirect, incidental, special or consequential damages arising out of the use of the Product, whether arising from an electrical issue or otherwise. Some States do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

5. Validity

This Limited Warranty is valid only in the USA.

6. Definitions

- a) "Normal Use" means ordinary consumer use under normal home conditions according to the instruction manual included with the Product. It does not include using the Product with a voltage converter.
- b) "Product" means **Opal**.
- c) "Warranty Period" means **2 years** from the date of purchase of the Product.

CLARISONIC PRODUCTS

Opal Sonic Infusion and our exclusive Anti-Aging Sea Serum are optimized to work in perfect synergy for best results. PBL is not responsible for any adverse results from using the Opal with any product other than the included Anti-Aging Sea Serum.

CLARISONIC SKIN CLEANSING

The Clarisonic Skin Cleansing systems gently cleanse the skin 6X better than by hand and prepare it to better absorb serums and moisturizers. Visit us at **clarisonic.com** for more information.





© 2012 Pacific Bioscience Laboratories, Inc.
All rights reserved. Clarisonic and Opal are
registered trademarks of Pacific Bioscience
Laboratories, Inc. (PBL).

Patents Pending

Pacific Bioscience Laboratories, Inc.
17275 Northeast 67th Court
Redmond, WA 98052 USA
Toll Free: 888-5-CLARISONIC or 888-525-2747

Every Opal Sonic Infusion System is hand
assembled in the USA.

PBL05530C